

AR reference images

Quest 3D DICOM 1.10.0 | Orthopractis

Print at actual size (100%)

Disable Fit to page. Each complete image below must measure 40 x 40 mm.


Keep the image flat. Do not crop, stretch, add marks or use a glossy cover.

Measure the scale bar with a ruler before using this sheet for placement.



RED REFERENCE

40 mm full-image width

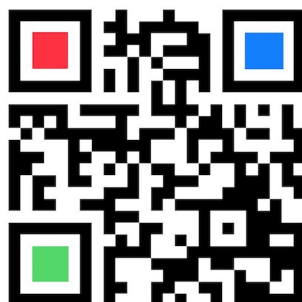

GREEN REFERENCE

40 mm full-image width



BLUE REFERENCE

40 mm full-image width



50 mm calibration bar

Use in the app

1. Choose the matching reference image in Target, then import content.
2. Leave the reference-image width at 40 mm in Tracking for this print.
3. Enter AR and show the complete image to the camera in good light.
4. Move slowly until the image is fully tracked and the origin appears.
5. Drag, pinch and twist to adjust content; save placement when aligned.
6. Keep the image visible for live tracking. Lock only a stationary placement.

The app recognises this exact artwork, not an arbitrary coloured sphere.

For another print size, enter its measured width in Tracking before restarting AR.